



動物用医薬品構造式一覧


PL動物薬 LC/MS Mix 1(Code:99056110)

PL動物薬 LC/MS Mix 2(Code:99056111)

TCs : テトラサイクリン系抗生物質

Sulfaguanidine


CAS : 57-67-0
MF : C7H10N4O2S
MW : 214.2
LogPow : -0.7

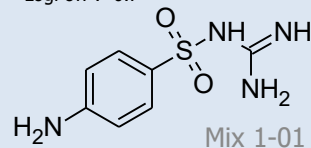

Mix 1-01

Sulfanilamide


CAS : 63-74-1
MF : C6H8N2O2S
MW : 172.2
LogPow : -0.6

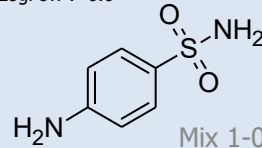

Mix 1-02

Sulfacetamide

CAS : 144-80-9
MF : C8H10N2O3S
MW : 214.2
LogPow : -1

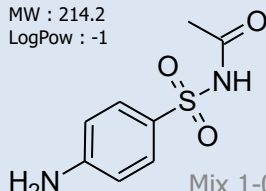

Mix 1-03

Sulfisomidine


CAS : 515-64-0
MF : C12H14N4O2S
MW : 278.3
LogPow : 1.2



Mix 1-04

Sulfadiazine


CAS : 68-35-9
MF : C10H10N4O2S
MW : 250.3
LogPow : -0.1

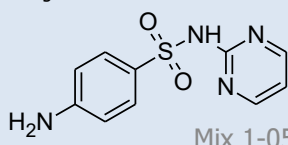

Mix 1-05

Sulfathiazole


CAS : 72-14-0
MF : C9H9N3O2S2
MW : 255.3
LogPow : 0.1

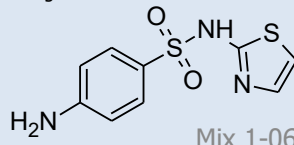

Mix 1-06

Sulfapyridine


CAS : 144-83-2
MF : C11H11N3O2S
MW : 249.3
LogPow : 0

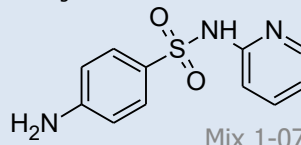

Mix 1-07

Diaveridine


CAS : 5355-16-8
MF : C13H16N4O2
MW : 260.3
LogPow : 1

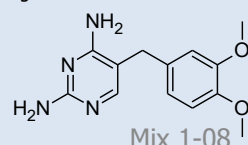

Mix 1-08

Sulfamerazine


CAS : 127-79-7
MF : C11H12N4O2S
MW : 264.3
LogPow : -1

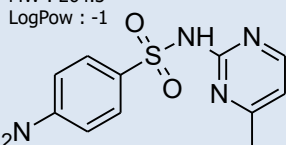

Mix 1-09

Trimethoprim


CAS : 738-70-5
MF : C14H18N4O3
MW : 290.3
LogPow : 0.9

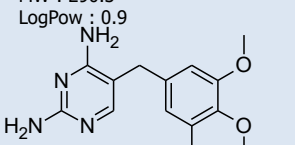

Mix 1-10

Sulfisozole


CAS : 73247-57-1(Na塩)
MF : C9H9N3O3S
MW : 239.3
LogPow : -

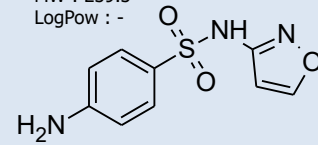

Mix 1-11

Ormetoprim


CAS : 6981-18-6
MF : C14H18N4O2
MW : 274.3
LogPow : 1.2

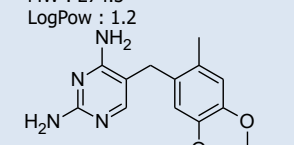

Mix 1-12

Sulfadimidine


CAS : 57-68-1
MF : C12H14N4O2S
MW : 278.3
LogPow : 0.3

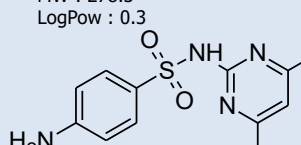

Mix 1-13

Sulfamonomethoxine


CAS : 1220-83-3
MF : C11H12N4O3S
MW : 280.3
LogPow : 0.8

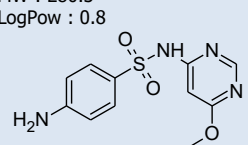

Mix 1-14

Sulfamethoxypyridiazine


CAS : 80-35-3
MF : C11H12N4O3S
MW : 280.3
LogPow : 0.3

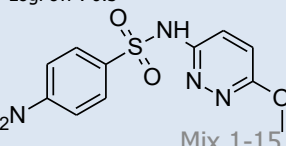

Mix 1-15

Sulfachlorpyridiazine


CAS : 80-32-0
MF : C10H9ClN4O2S
MW : 284.7
LogPow : 1

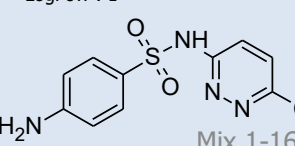

Mix 1-16

Sulfamethoxazole


CAS : 723-46-6
MF : C10H11N3O3S
MW : 253.3
LogPow : 0.9

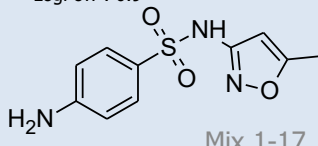

Mix 1-17

Sulfadoxine


CAS : 2447-57-6
MF : C12H14N4O4S
MW : 310.3
LogPow : 0.7

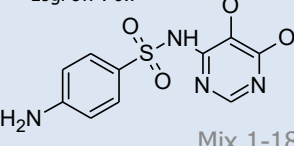

Mix 1-18

Sulfatroxazole


CAS : 23256-23-7
MF : C11H13N3O3S
MW : 267.3
LogPow : 1.3

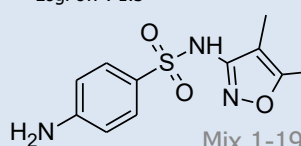

Mix 1-19

Sulfaethoxypyridazine


CAS : 963-14-4
MF : C12H14N4O3S
MW : 294.3
LogPow : 0.7

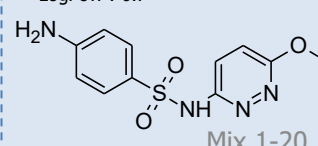

Mix 1-20

Sulfisoxazole


CAS : 127-69-5
MF : C11H13N3O3S
MW : 267.3
LogPow : 1

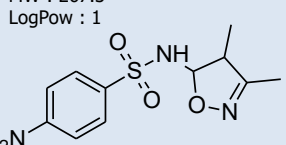

Mix 1-21

Pyrimethamine


CAS : 58-14-0
MF : C12H13ClN4
MW : 248.7
LogPow : 2.7

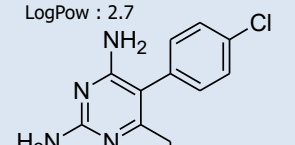

Mix 1-22

Sulfabenzamide


CAS : 127-71-9
MF : C13H12N2O3S
MW : 276.3
LogPow : 0.7

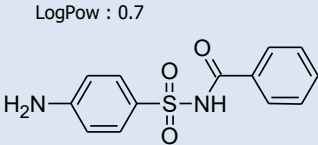

Mix 1-23

Sulfadimethoxine


CAS : 122-11-2
MF : C12H14N4O4S
MW : 310.3
LogPow : 1.6

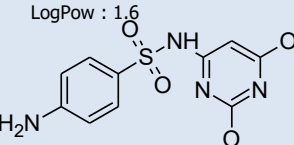

Mix 1-24

Sulfaquinolaxine


CAS : 59-40-5
MF : C14H12N4O2S
MW : 300.3
LogPow : 1.7

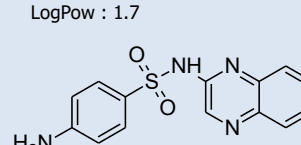

Mix 1-25

Sulfabromomethazine


CAS : 116-45-0
MF : C12H13BrN4O2S
MW : 357.2
LogPow : 1

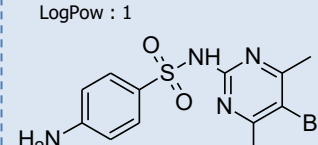

Mix 1-26

Sulfantran


CAS : 122-16-7
MF : C14H13N3O5S
MW : 335.3
LogPow : 2

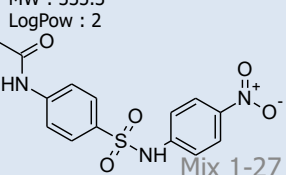

Mix 1-27

Marbofloxacin


CAS : 115550-35-1
MF : C17H19FN4O4
MW : 362.4
LogPow : -0.5

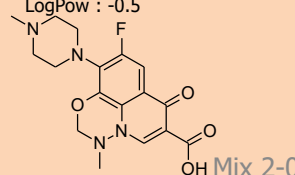

Mix 2-01

Norfloxacin


CAS : 70458-96-7
MF : C16H18FN3O3
MW : 319.3
LogPow : -1

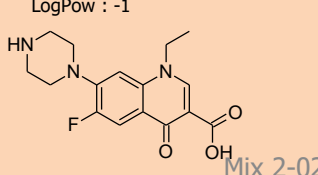

Mix 2-02

Ofloxacin


CAS : 82419-36-1
MF : C18H20FN3O4
MW : 361.4
LogPow : -0.4

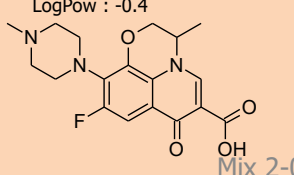

Mix 2-03

Ciprofloxacin


CAS : 85721-33-1
MF : C17H18FN3O3
MW : 331.3
LogPow : -1.1

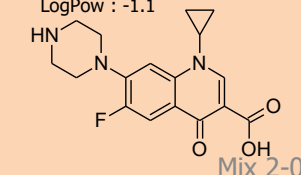

Mix 2-04

Danofloxacin


CAS : 112398-08-0
MF : C19H20FN3O3
MW : 357.4
LogPow : -0.3

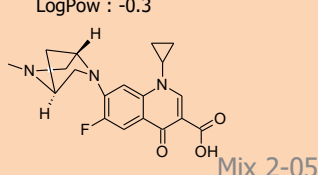

Mix 2-05

Enrofloxacin


CAS : 93106-60-6
MF : C19H22FN3O3
MW : 359.4
LogPow : -0.3

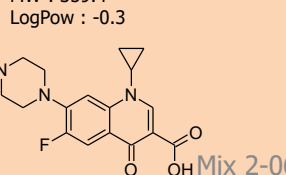

Mix 2-06

Orbifloxacin


CAS : 113617-63-3
MF : C19H20F3N3O3
MW : 395.4
LogPow : 0.9

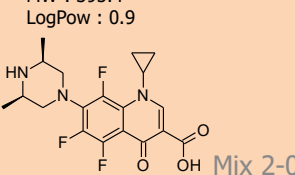

Mix 2-07

Sarafloxacin


CAS : 98105-99-8
MF : C20H17F2N3O3
MW : 385.4
LogPow : 0.3

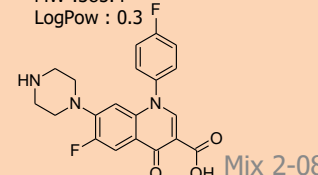

Mix 2-08

Difloxacin


CAS : 98106-17-3
MF : C21H19F2N3O3
MW : 399.4
LogPow : 1.6

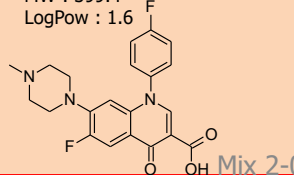

Mix 2-09

Miloxacin


CAS : 37065-29-5
MF : C12H9NO6
MW : 263.2
LogPow : 1.8

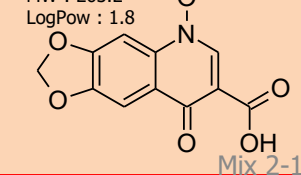

Mix 2-10

Oxolinic acid


CAS : 14698-29-4
MF : C13H11NO5
MW : 261.2
LogPow : -0.2

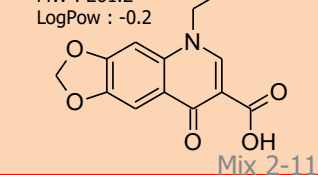

Mix 2-11

Nalidixic acid


CAS : 389-08-2
MF : C12H12N2O3
MW : 232.2
LogPow : 1.4

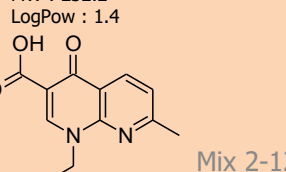

Mix 2-12

Flumequine


CAS : 42835-25-6
MF : C14H12FN3O3
MW : 261.3
LogPow : 2.9

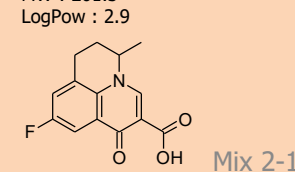

Mix 2-13

Piromidic acid


CAS : 19562-30-2
MF : C14H16N4O3
MW : 288.3
LogPow : 1.3

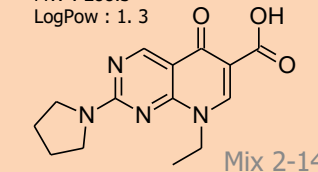

Mix 2-14

Tetracycline


CAS : 60-54-8
MF : C22H24N2O8
MW : 444.4
LogPow : -2

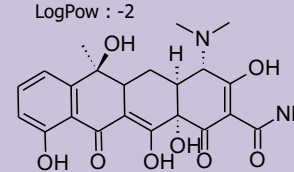

TCs-01

Oxytetracycline


CAS : 79-57-2
MF : C22H24N2O9
MW : 460.4
LogPow : -1.6

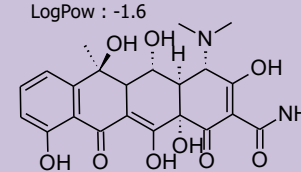

TCs-02

Chlortetracycline

CAS : 57-62-5
MF : C22H23ClN2O8
MW : 478.9
LogPow : -1.3

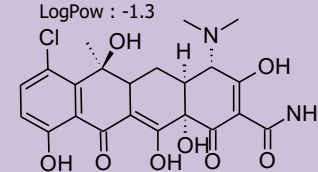

TCs-03

株式会社 アイスティサイエンス

<http://www.aisti.co.jp>

Tel:073-475-0033 Fax:073-497-5011

E-mail:as@aisti.co.jp